


NEO₂ BLEND: Air / Oxygen Blender for the NICU

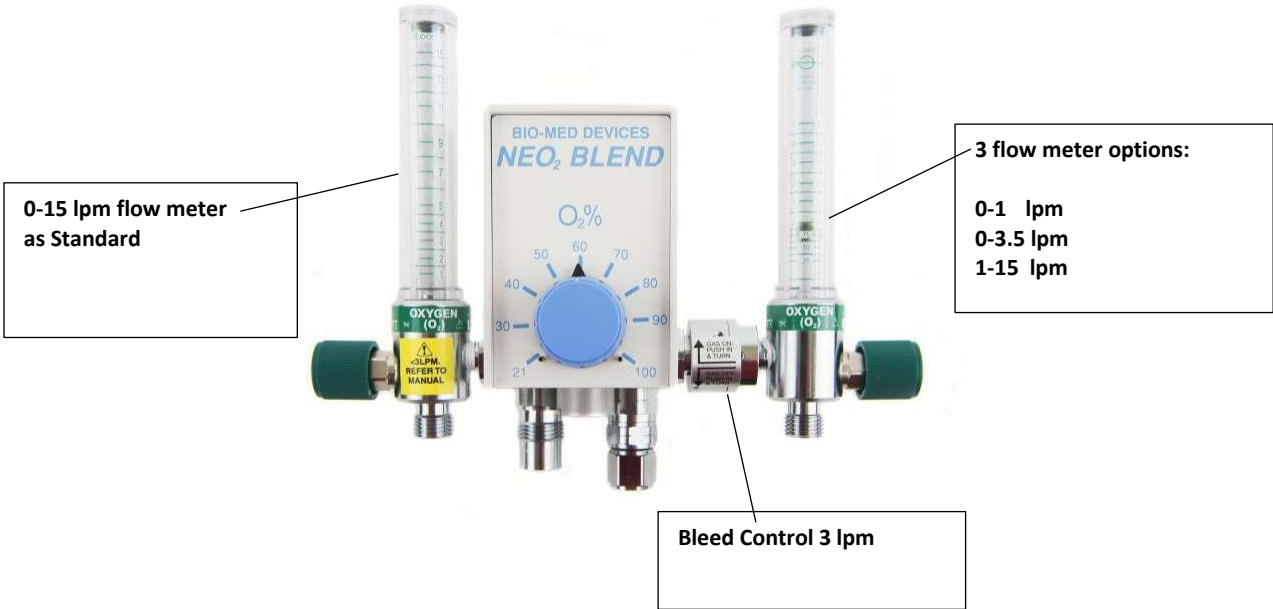

Ideal model for NICU Application

Key features

Comes with 0-15 lpm flow meter close mounted to the left-hand side

This makes a perfect spot to attach a resuscitation bag and keeps the footprint as small as possible and minimizes the amount of wall space required position and use as desired

3 Different options configured with the Right-side flow meter depending on the needs and practice of the unit

Right side comes installed with either 1) 0-1 lpm 2) 0-3.5 lpm or 3) 1-15 lpm flowmeter depending on the model you choose, which allows for use with a nasal cannula, CPAP, etc.

Unique feature on the right-side BLEED Control (3 lpm)

The right-side flow meter is installed on an integrated bleed control switch: When the right side is not needed, the bleed flow (3 lpm) can be disabled by pushing and rotating the flow meter forward when needed, simply return it to the upright position and use as desired

Standard input fittings are DISS for Air and Oxygen

But, NIST fittings on the blender is also available.

NEO2 BLEND: Air / Oxygen Blender for the NICU

Technical Specifications

Oxygen Concentration Range	21 to 100 %
Number of ports	2 or 3 (2 occupied by flow meters, bottom port present on 2003FFA model)
Left Flowmeter Range	3-15 lpm (1-15 if bleed is on)
Right Flowmeter Range	0-1.0 lpm 0-3.5 lpm 1-15 lpm (with 3 lpm bleed on- flowmeter in upright position)
Accuracy	+/- 3% full scale
Bleed flow	3 lpm (when activated)

Ordering information

2003FF1	NEO2 Blend with	15 lpm	1.0 lpm	O2 Green
2003FF	NEO2 Blend with	15 lpm	3.5 lpm	O2 Green
2003FF15	NEO2 Blend with	15 lpm	15 lpm	O2 Green
2003FFA	NEO2 Blend with	15 lpm	3.5 lpm	O2 White
2003FFEU	NEO2 Blend with	15 lpm	3.5 lpm	O2 White
2003FFE5	NEO2 Blend with	15 lpm	15 lpm	O2 White
2003FFEA	NEO2 Blend with	15 lpm	3.5 lpm	O2 White and bottom output

Dimension & Weight

Height	Width	Depth	Weight
6.64" (16.9 cm)	7.98" (20.3 cm)	4.56" (11.6 cm) w/ flow meter on 6.69" (17 cm) w/ flow meter off	3.17 lbs. (1.44 kg)

Optional Accessories

High Pressure Hose Kit	2', 3' or 10' long
Bracket	Wall Mount, Pole Mount, Rail, headwall configuration
Quick Connect Fittings	Ohmeda, Chemtron, Quick Connect fittings

**This data is subject to change without a prior notice.*

